

haakaa®

SILICONE SPILLGUARD™ INSERT

Product Manual

MHK243 E10319

Product Name: Silicone SpillGUARD™ Insert**Material:** Food-Grade Silicone**Dimensions:** 49 x 42 x 52mm

Product Description:

Prevent accidental spills without interrupting your pumping sessions using Haakaa's Silicone SpillGUARD™ Insert! It sits snugly in the necks of all Haakaa Gen. 1, 2, 2* & 3 pumps while you're pumping - without affecting suction. Compact, durable, and portable, it even includes an easy-grip tab for simple insertion and removal. Because this innovative valve stopper stays in the pump while you're using it, you'll never have to hunt around to find a way to keep your milk secure after you've finished again! Made from 100% food-grade silicone, this stopper is the ideal way to give yourself peace-of-mind while expressing, preventing tears over spilled milk!


Features:

- Fits snugly into the neck of any Haakaa Gen. 1, 2, 2* or 3 pump without affecting suction.
- Creates a secure seal to prevent accidental spills during pumping.
- Non-slip tab for simple insertion, alignment and removal.
- Compact, portable & durable to help you pump securely anytime, anywhere.
- Safe in microwaves, boiling water, dishwashers and sterilisers.
- BPA, PVC and phthalate-free.


Cleaning and Care:

Remove stopper from pump before cleaning. Wash before and after each use. Although dishwasher safe, we recommend hand washing in warm soapy water and rinsing and drying thoroughly. Only soft bristle brushes or soft sponges should be used to clean this product. Hard scourers may scratch the surface. Do not use any bleach-based cleaners or tablets to sterilise or clean your Haakaa products. To sterilise, use a steam steriliser or boil in water for 2-3 minutes. Do not use UV sterilisers to sterilise this product. Store product in a cool dry place and avoid direct sunlight.


Instructions:

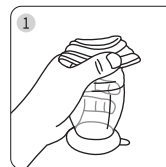


Wash your hands and all components before use. Fold down the flange at the narrowest part of the neck.

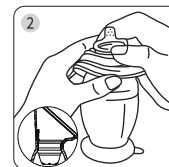

Face the non-slip tab upwards and align with the spillway of the pump. Firmly insert the stopper into the neck of the pump, ensuring there are no gaps.

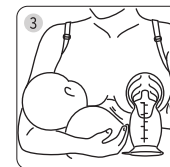

Fold flange back and use as usual once stopper is inserted. Squeeze the base, centre the opening over your nipple, then release.

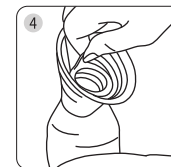

To remove, simply lift the tab and pull up.

PLEASE NOTE: The stopper is designed with small holes to stabilise pressure. To prevent leaks, avoid tipping the pump upside down, even with the stopper inserted.

Caution:

- Check product condition regularly. Replace if product shows any signs of damage.
- Do not store near any sharp objects.
- This product is not a toy. Do not use for anything other than its intended purpose.
- Keep out of the reach of children and pets.
- As milk quantity increases in the pump, additional squeezing of the pump body may result in air or milk spurting out. Remove stopper, transfer expressed milk to a clean, sterilised container then reinsert stopper to continue pumping without leaks.

PLEASE NOTE: This product is an accessory for our Silicone Breast Pumps and will not work by itself. Always follow your country's guidelines and recommendations for safely expressing and storing breastmilk.