尺寸:
210x200mm

工艺要求:
70g双胶纸

半成品号:
E10945

haakaa®

LATCHFLOW SYSTEM™ NIPPLE SHIELD
Product Manual English

MHK265 E10945

Product Name: LatchFlow System™ Nipple Shield

Materials:
Nipple Shield – Medical-Grade Silicone
Storage Frame & Case – PP
Cleaning Brush – Stainless Steel, Nylon, PP

Dimensions:
Nipple Shield (S) – 131 x 63 x 77 mm
Nipple Shield (M) – 131 x 63 x 80 mm
Nipple Shield (L) – 131 x 63 x 81 mm
Storage Frame – 79 x 69 x 49 mm
Cleaning Brush – 16 x 9 x 150 mm
Storage Case – 117 x 81 x 77 mm

Opening Size:
S – 17 mm
M – 20 mm
L – 24 mm
Connection Opening: 3.8 mm

Nipple Shield Standard: EN 14350

Intended Use:

This nipple shield is intended to provide supplemental milk to babies while they feed at the breast when connected to a compatible milk source. Use only as recommended by a lactation professional.

Product Description:

Designed to support breastfeeding through at-breast supplementation, this double-layered nipple shield helps deliver additional nutrition while baby feeds naturally at the breast. Built-in channels helps guide supplemental milk toward the nipple.

The ultra-thin, cut-out design enhances skin-to-skin contact while offering gentle protection against irritation from baby's sucking and biting. Soft, flexible medical-grade silicone ensures a comfortable fit, while small openings allow mum's let-down to pass through for a more intuitive feeding experience. Finally, keep your shield tidy and organised with the included storage case and shape-retaining frame!

Features:


- Provides extra supplemental milk to baby while allowing mum and baby to continue breastfeeding.
- Helps mothers continue breastfeeding and maintain skin-to-skin contact while lactation is being established.
- Fits Haakaa LatchFlow System™ pouches, dispensers, supplemental feeding systems, or oral feeding syringes (none included) via the 3.8 mm opening.
- Double-layered silicone shield with built-in channels helps guide supplemental milk toward the nipple.
- Ultra-thin shield has unique cut-outs to allow closer contact for a more natural feeding experience.
- Made of skin-friendly 100% medical-grade silicone.
- Openings on the nipple shield allow milk to pass through from mum's nipple.
- Includes a shape-retaining frame for the nipple shield and a storage case to keep all parts neatly organised.
- Safe in dishwashers, sterilisers and boiling water with a dedicated brush for thorough cleaning.
- BPA and phthalate-free.

Cleaning and Care:

◎ Before first use, disassemble and clean all parts of this product before placing silicone components in boiling water for 5 minutes to sterilise. Boil plastic part (storage frame & case) for 2-3 minutes to sterilise (do not boil any longer as this may impact the lifespan of plastic parts). Wash as soon as possible after each use to prevent milk residue drying.
◎ To clean the nipple shield and tubing, press the centre dome through from the other direction to separate it from the outer flange. Flush with warm, soapy water and use the included cleaning brush to gently scrub all shield surfaces, both inside and out. Gently insert the cleaning brush into the tubing and move it back and forth to clean the inside. Rinse then leave inverted to thoroughly dry.
◎ To sterilise all parts except the cleaning brush, boil in water for 2-3 minutes. Brush **NOT** suitable for microwaves, sterilisers or boiling water due to stainless steel components. Do not use any bleach-based agents or sterilising tablets to clean this product. Do not use a UV steriliser to clean this product as it may impact its lifespan.
◎ Only soft bristle brushes or soft sponges should be used to clean this product, as hard scourers may scratch the surface. Store in a cool, dry place and avoid direct sunlight.

PLEASE NOTE: Always follow your country's guidelines and recommendations for safely expressing and storing breastmilk.

Find Your Perfect Fit:

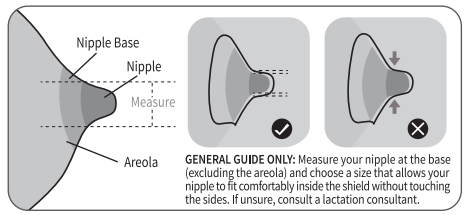


GENERAL GUIDE ONLY: Measure your nipple at the base (excluding the areola) and choose a size that allows your nipple to fit comfortably inside the shield without touching the sides. If unsure, consult a lactation consultant.

Shield Opening Size	Recommended Nipple Size
S: 17 mm	13-16 mm
M: 20 mm	16-19 mm
L: 24 mm	19-23 mm


NOTE: If between sizes, start with the smaller shield. If pinching, rubbing or discomfort occurs, consider selecting the next size up. If there is too much movement inside the shield, or your nipple does not sit comfortably within the opening, try the next size down. You may find it helpful to measure your nipple again during use, as nipple size and shape can change during breastfeeding. We also recommend remeasuring at 6-8 weeks postpartum, or any time feeding comfort changes. For further information, contact your local IBCLC, lactation consultant or maternal health provider.

Instructions:


Ensure each part is correctly fitted and securely connected before use.

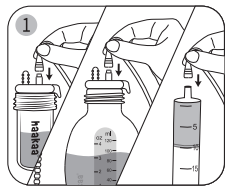


Wash all parts and your hands before use. Fill dispenser, pouch or feeding syringe with supplemental milk - none included. Push the silicone cap and connector firmly into the pouch opening.

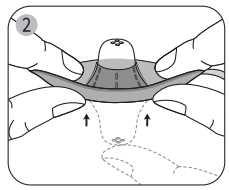


Ensure the two layers of the shield overlap correctly. Carefully turn brim upwards.

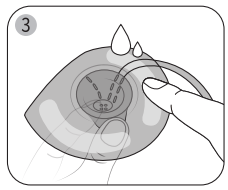


Place a few drops of drinking-quality water or breastmilk on the underside of the brim to help the shield stick to your skin.

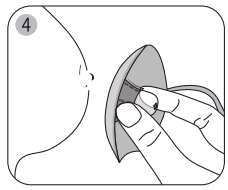

Centre the shield over your nipple, enclosing as much of your nipple as possible.

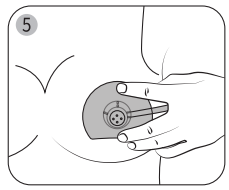

Press the shield down onto your skin to turn the brim back down. The base should sit flush against your skin.

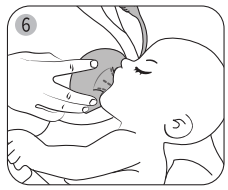

With the shield in place, gently introduce your baby to your nipple. Position baby on the shield to allow milk from the container to flow into their mouth.

NOTE: Dispenser, pouch and oral feeding syringe not included. For supplemental feeding, connect the nipple shield to a compatible haakaa LatchFlow System™ dispenser or pouch, or to a compatible oral feeding syringe. Ensure the connector fits securely into the 3.8 mm opening before use.

Caution:

For your child's safety and health
WARNING!

◎ Continuous and prolonged sucking of fluids will cause tooth decay.
◎ Always check food temperature before feeding.
◎ Throw away at the first signs of damage or weakness.
◎ Keep components not in use out of the reach of children.
◎ Never attach to cords, ribbons, laces or loose parts of clothing. The child can be strangled.
◎ Always use this product with adult supervision. Only adults should attempt to feed babies from this product.
◎ Clean and dry all tubing thoroughly after each use to prevent milk residue buildup.
◎ This is a single-user product. Use by more than one user without adequate cleaning may present a health risk and should be avoided.
◎ We recommend replacing feeding shield regularly.
◎ This product should only be used under the advice of a lactation professional. Always consult a medical professional for further advice.
◎ This product is not a toy. Do not use this product for anything other than its intended use.
◎ Check all tubing, connectors and parts are securely assembled before each use.
◎ Never use feeding teats as a soother.
◎ Do not leave a feeding teat in direct sunlight or heat, or leave in disinfectant ("sterilising solution") for longer than recommended, as this may weaken the teat.
◎ If using third-party feeding components, ONLY use ones that securely fit the 3.8 mm opening. Do not force incompatible parts into the opening.
◎ Do not forcefully squeeze or shake at any time.